


Libris Mobile Fall Detector

DISCLAIMER: MHM will do whatever we can to support this installed equipment. We take responsibility for things in our control and will support you and the installation within our normal operating procedures. We cannot assume responsibility for monitoring centers, EMS, equipment misuse, weather disturbances, phone networks of any kind, cell phone towers or coverage, or other events or processes beyond our control.

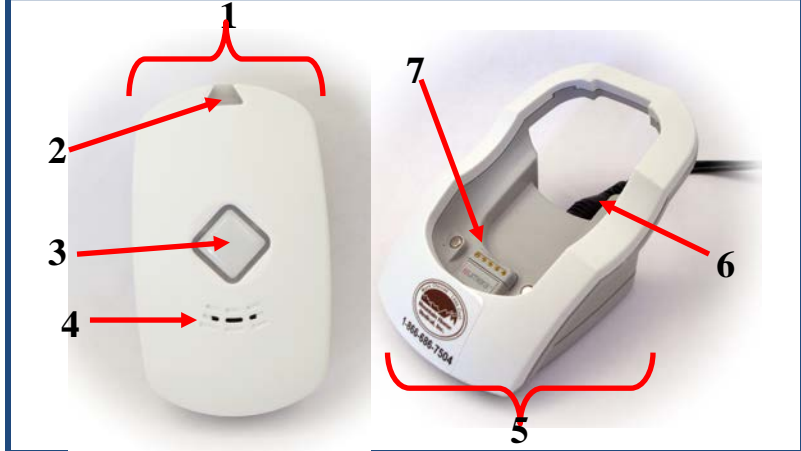

<u>Parts of the Libris</u>	<u>Part Names / Descriptions</u>
	<ol style="list-style-type: none"> 1. Mobile Device 2. Cellular Signal Strength 3. Battery Status / Call Button 4. Microphone 5. Charging Cradle 6. Power Cord 7. Charging Contacts

Using the Libris







1. **If you are doing a self- install and were not visited by a Mountain Home Medical Technician**, call MHM at 1-866-686-7504 and inform them you have received your Libris device and would like to activate it. If you don't call the office for activation BEFORE turning it on and testing, it will not function correctly.
2. Remove your Libris from the box. Locate the charging cradle and connect the small round end of the power cord into the back of the charging cradle. Plug the other end of the cord into a working electrical outlet that is not controlled by a light switch. It is recommended to place the charging cradle near your bed so you can access it at night if needed.
3. **Charge the unit.** Place the mobile device into the charging cradle bottom first. The device will power on by itself and the indicator lights will flash different colors. It will say "We are activating your personal mobile device. We will notify you when activation is complete". This registration process can take up to 30 minutes. When complete, it will say "Activation is now complete. We are resetting your personal mobile responder so it is ready for your use. One moment please". The unit should be charged at least 4 hours before use. Failure to fully charge the Libris will significantly shorten its battery life. **If the battery dies, the Libris will not function until recharged and powered back on.**
4. **Test the Libris.** Press and hold down the Call Button in the center of the mobile device until it lights up (blue color). You will hear a voice prompt that it is calling your response team. Allow the device to connect to the response center and inform them you are testing your new device.

A NOTE ABOUT GPS: GPS technology in a device this small is not 100% accurate. Buildings, mountains, weather, and electrical interference can all effect how it operates. This unit can detect the general vicinity of where you are, but sometimes may be off by a few houses or a street. It is best to always be aware of your surroundings so if you need help, you can tell the operator where you are.

Indicator Lights on the Libris Mobile Device

Light	Appearance	Status
 Battery Status	Green	Good Battery Life – no charging required
	Amber	Moderate Battery; charge within 2 hours
	Red	Low Battery; requires charging now
	Color Cycling	Unit has turned on and is charging (unit was turned off)
	Off (gray)	Unit is in standby mode or turned off
	Blue	Unit has been activated and is calling the response center
 Cellular Signal	Green	Strong ATT signal (5 bars)
	Amber	Moderate ATT signal (3 bars)
	Red	No ATT signal, unit will not function
	Off (gray)	Unit is in standby mode or turned off

Features of the Libris

 Extended battery life (up to 36 hours in areas with strong ATT cell service)	 Water / Splash Resistant (Do not submerge or wear while swimming)
 Easy to wear (includes neck lanyard, and belt or purse clip)	 Smart Fall Detection (unit learns patient's patterns and gait to help prevent false alarms)
 24/7/365 Response Center (Trained staff ready to take your call at any time)	 Voice Prompts (instead of beeps, the device has voice prompts to tell you its status and what it is doing)

Basic Operation

Turning Libris On	Place Libris in the charging cradle and it will turn on automatically. If the cradle is not available, press and hold the Call Button until the light flashes. It will play music while it starts up. Once Connected, Libris will say " <u>Your personal mobile responder is ready</u> ".
Turning Libris Off	Press the Call Button twice, holding down the second press for 7 seconds until Libris plays a message that you have requested to turn it off. Press the call button once more to confirm. If done correctly, Libris will say " <u>Powering off now. Goodbye.</u> "
Charging the Libris	Remove any clips or lanyards and place the Libris bottom first into the charging cradle. Make sure the charging pins on the cradle line up with the metal dots on the back of the unit. The Libris must be removed from the charging cradle daily to allow the battery to drain before it can be re-charged. If left on the charger for extended periods, the battery will no longer take a charge and the unit must undergo a factory reset in order to charge again.
Using Fall Detection	The technology used in the Libris to detect falls is dependent on several sensors, accelerometers, and computer programs. The libris should be worn at all times for the first 3 weeks of use so the program can 'learn' how you walk, sit down, stumble, lay down, etc. so it will have fewer false alarms later. During this break-in period, the Libris may have some false alarms, and it may not activate for some falls because it is still learning. No fall detection technology is 100% accurate; if you fall and need help, don't wait for the unit to activate – press the call button.

TECHNICAL SUPPORT: If you have any questions, concerns, or your unit stops working, please call our office at 1-866-686-7504 and our friendly office staff would be happy to help you.

-----TO BE LEFT WITH CLIENT-----